


15" Compact Freezer

Model # MS15FAS6RW

- Front vented, permitting it to be built into casework or installed under lab benches without the need for additional ventilation clearance
- Mechanical thermostat factory preset at -20 °C
- Accurate temperature control of contents within +3°/
 -3°C of set point.
- Recommended ambient temperature 18.3°C to 32.2°C (65°F to 90°F)
- Non-CFC 245FA insulation in cabinet and door seals in the cold
- 22 gauge textured steel cabinet and door exterior
- Forced air condenser system efficiently removes heat and allows this freezer to be used on a counter top or built into casework without the loss of cooling performance
- Powerful compressor for fast pull-down time and temperature recovery after door opening or commodity loading
- "Box Style" door for easy mounting of inventory management system hardware
- Discrete top edge mounted door handle
- Cam-assisted door close hinge system to reduce "door-ajar" events
- Durable white powder coat paint finish
- CSA Listed and warrantied for commercial use

Rough In Dimensions

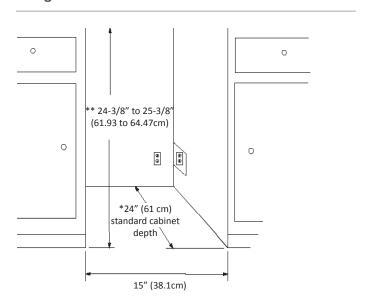



Ordering Details	
Solid White Cabinet and Door with Right Hinge	MS15FAS6RW

Width	14 7/8"	
Interior Finish	White	
Total Capacity	1.5 Cubic Feet; 42.5 Liters	
Toe Grill Finish	Black	
Hinge	Nylon bushing with hinge pin	
Handle	Discrete	
Lock	No	
Temperature Range	-4°F (-20°C)	
Electrical Requirements	115V / 60Hz / 3.3A	
Length of Power Cord	5'	
Product Dimensions	14 ⁷ / ₈ " W x 24 ¹ / ₄ " H x 22 ¹ / ₈ " D	
Shipping Weight	90 lbs; 122 lbs with skid	
Finishes	White	

Product Dimensions

