

# 24" General Purpose All Freezer

Model # MS24FAS2RW



# Door-front mounted MicroSentry<sup>™</sup> II Scientific refrigeration controller provides:


- Accurate electronic temperature control of contents within +2°/ -2°C of set point.
- Recommended ambient temperature 18.3°C to 32.2°C (65°F to 90°F)
- Large, easy-to-read digital temperature display with resolution to 0.1°C
- User selectable temperature set point between -12°C and -25°C, factory preset to -20°C
- User selectable high and low temperature alarm values
- User selectable alarm delay time period
- Visual and audible alarms for over temperature, under temperature, and door ajar events
- Alarm silence and automatic ring back functions
- Keypad lock-out to prevent unintended control inputs
- Basic high and low temperature logging with
- reset capability
- Visual diagnostic error codes



#### **Additional Features**

- Outstanding cooling capacity and temperature stability
- One full depth wire shelf
- Front vented through the toe space grill for true built in installations
- Powerful compressor for fast pull-down time and temperature recovery after door opening or content loading
- Manual defrost eliminates temperature excursions and dehydration
- "Box Style" door for easy mounting of inventory management system hardware
- Discrete top edge mounted door handle
- Cam-assisted door close hinge system to reduce "door-ajar" events
- New triple-seal energy efficient door gasket system seals in the cold, reducing energy consumption
- 1 diameter probe port in rear cabinet wall (standard)
- 4 leveling legs for easy height adjustment
- Extra-long, jacketed commercial power cord with hospital grade plug to meet stringent hospital and clinical use standards
- Durable white powder coat paint finish
- CSA listed and warrantied for commercial use

## **Rough In Dimensions**

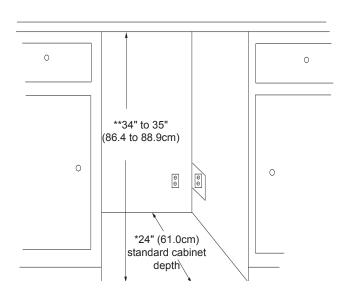


## **Finish/Door Options**



<b>Ordering Details</b>	
Solid White Cabinet and Door with Lock, Right Hinge	MS24FAS2RW

#### **Product Dimensions**

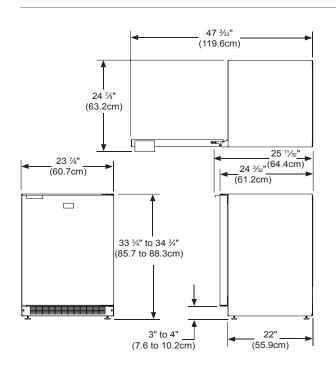



Width	23 7/8"
Interior Finish	White
Total Capacity	4.5 Cubic Feet; 127 Liters
Toe Grill Finish	Black
Hinge	Nylon bushing with hinge pin
Handle	Discrete
Lock	Yes

Temperature Range	10°F to -13°F (-12°C to -25°C)
Electrical Requirements	115V / 60Hz / 3.3A
Length of Power Cord	5'
Product Dimensions	23 <sup>7</sup> /8" W x 34" H x 24" D
Shipping Weight	120 lbs; 130 lbs with skid
Finishes	White
Agency	CSA Listed