

HERMLE


Benchmark

Z306

Universal Centrifuge



- Universal, 9 available rotor options
- Compact, less than 14 inches wide
- Advanced programming with auto rotor recognition
- Swing out rotor for 8 x 50ml or 12 x 15ml

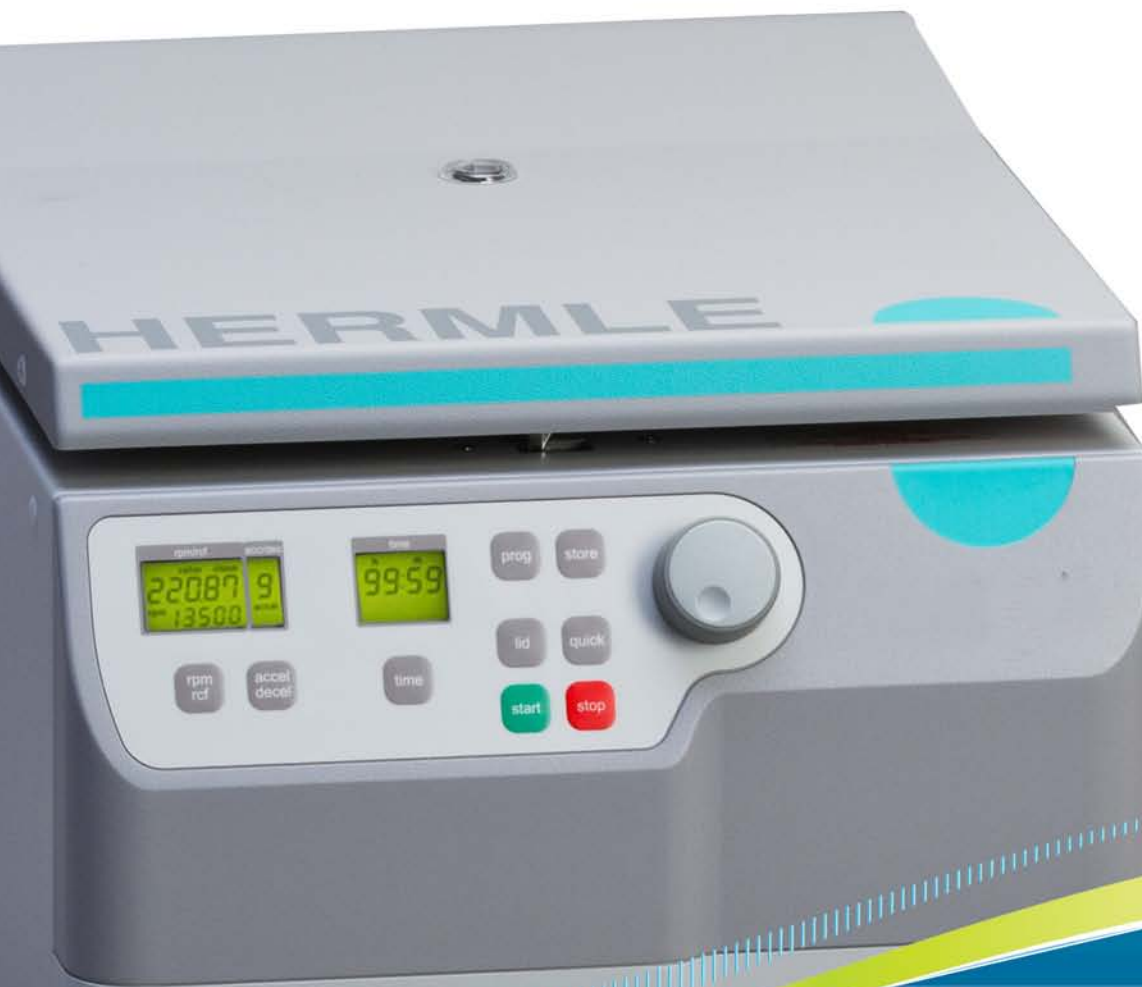

Z306 Universal Centrifuge

Z306 Summary

Speed up to 14,000 rpm
G-Force up to 18,624 x g
Swing Out Capacity:
8 x 50 ml, 12 x 15 ml,
6 x microplates

Angle Rotor Capacity
30 x 1.5/2.0 ml
6 x 50 ml
30 x 15 ml

Z306 Features

Microprocessor Control
Over-speed protection
Automatic rotor identification
Program Storage (99)
Momentary spin button
Imbalance detection
Adjustable Accel/Decel
Audible end of run alarm

Maintenance free induction drive
Electrical safety lid-lock
Quiet operation, less than 63 DB
Splash proof display

Z306 Specifications

Speed Range: 200 to 14,000 rpm
Max. RCF: 18,624 x g
Max Capacity: 4 x 100 ml
Timer: up to 99 hr, 59 min.
Dimensions: 13.9 x 18.6 x 13 in.
(W x D x H) 35.5 x 47.4 x 33 cm.
Weight: 66 lbs (30 kg)
Electrical: 120 or 230V, 50/60Hz
Power Draw: 240W



The Z306 is a compact and economical centrifuge with a series of rotor options making it well suited for clinical, industrial, and research laboratories.

Despite the impressive capacity of the 4 x 100ml rotor (capable of supporting 8 x 50ml or 12 x 15ml with the use of rotor inserts), the Z306 is surprisingly compact with a footprint of just over 14 x 18 inches. An optional rotor is also available for microplates (2 x 3 microplates), positioning this model as one of the smallest, and least expensive centrifuges capable of high speed microplate centrifugation.

Advanced microprocessor controls allow for instant rotor identification, which automatically converts the rotational speed to the relevant g-force (rcf), preventing the user from setting a speed above the capabilities of the installed rotor. In addition to retrieving saved programs, the easy to use front panel can adjust accel/decel rates, and change the end of run alarm tone.

Ordering Information:

Z306* Z306 Centrifuge (without rotor), 120V
*115V with US Plug. For 230V (EU Plug), please order -E.

Universal Swing Out Rotors:



Item: Z306-100
4 x 100ml Swing Out Rotor (90°)
Max Speed: 4,500 rpm / 3,350 x g
Accel: 14 sec. / Decel: 19 sec.

Inserts:	Tube	qty/ insert	qty./ rotor
Z306-100-A100:	100ml	1	4
Z306-100-A85:	85ml	1	4
Z306-100-A50:	50ml conical	1	4
Z306-100-A15:	15ml conical	2	8
Z306-100-A10:	10ml	7	28
Z306-100-A57:	5/7ml	7	28
Z306-100-A20:	1.5/2.0ml	10	40

*Inserts sold in 2/pk.
Rotor accepts 4 inserts (2 x 2/pk.)



Item: Z326-100H
4 x 100ml Swing Out Rotor with hermetically sealed buckets (90°)
Max Speed: 4,000 rpm / 2,486 x g
Accel: 12 sec. / Decel: 12 sec.

Inserts	Tube	qty/ insert	qty./ rotor
Z326-100H-A100:	100ml	1	4
Z326-100H-A85:	85ml	1	4
Z326-100H-A50:	50ml conical	1	4
Z326-100H-A15:	15ml conical	2	8
Z326-100H-A10:	10ml	4	16
Z326-100H-A57:	5/7ml	4	16
Z326-100H-A20:	1.5/2.0ml	5	20

*Inserts for Z326-100H-AC Rotor, sold in 2/pk.
Rotor accepts 4 inserts (2 x 2/pk.)

Item: Z326-100H-AC

Same as above,
without buckets

Buckets	Tube	qty/ Bucket	qty./ rotor
Z326-100H-B50	50ml	2	8
Z326-100H-B15	15ml	3	12
Z326-100H-BUK	Hermetically Sealed Buckets (for use with Z326-100H-XXX inserts above)		

*Buckets for Z326-100H-AC rotor sold in 2/pk.
Rotor accepts 4 buckets (2 x 2/pk.)



Z326-100H-AC



B50

B15

Additional Rotors (Low Speed):



Item: Z326-0650-LS
6 x 50ml Angle Rotor (40°)
Max Speed: 6,000 rpm / 4,427 x g
Accel: 11 sec. / Decel: 14 sec.

Inserts/Adapters

Z206-50-A15: 15ml (conical or round)
Z206-50-A57: 5-7ml (12-13mm)
Z206-50-A57C: Cushion for 5ml tube



Item: Z326-02MP
Swing Out Rotor 2 x 3 microplates (90°)
Max Speed: 4,500 rpm / 2,716 x g
Accel: 23 sec. / Decel: 20 sec.

Inserts

BT1500-A1520: 20 x 1.5ml
BT1500-A05: 20 x 0.5ml
Z326-MPCA: 3 x chamber slide and dish carrier 2/pk

New optional accessory!!!

Z326-MPCA

3 x chamber slide and dish carrier



Chamber Slides

35mm Dishes

60mm Dishes

Fits all Hermle microplate carriers!!!



Item: Z326-1215
12 x 15ml Angle Rotor (32°)
Max Speed: 6,000 rpm / 4,427 x g
Accel: 11 sec. / Decel: 14 sec.
Inserts/Adapters
C3100-ADP: 5-7ml (12-13mm)



Item: Z326-3015
30 x 15ml Angle Rotor (35°)
Max Speed: 4,500 rpm / 2,830 x g
Accel: 19 sec. / Decel: 20 sec.
Inserts/Adapters
C3100-ADP: 5-7ml (12-13mm)

Additional Rotors (High Speed):



Item: Z326-2420
24 x 1.5/2.0ml Angle Rotor (40°)
Max Speed: 14,000 rpm / 18,624 x g (Z326)
Accel: 25 sec. / Decel: 36 sec.



Item: Z326-2420H
24 x 1.5/2.0ml Angle Rotor (40°) Hermetically Sealed Lid
Max Speed: 14,000 rpm / 18,624 x g
Accel: 25 sec. / Decel: 36 sec.



Item: Z326-3020H
30 x 1.5/2.0ml Angle Rotor (40°) Hermetically Sealed Lid
Max Speed: 12,000 rpm / 15,131 xg
Accel: 27 sec. / Decel: 24 sec.

Also Available from Benchmark:

Z446 Series High Capacity Centrifuges

Swing Out Rotor Accepts up to:

- 56 x 15ml conical tubes
- 28 x 50ml conical tubes
- 12 x microplates

